



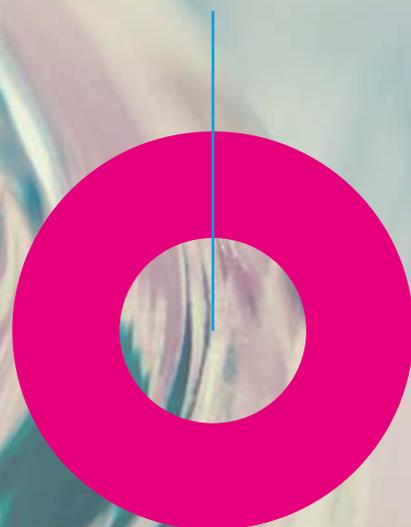

Your needs, our innovations

ΔΙΑ•Μ
современная лаборатория

www.dia-m.ru
заказ on-line

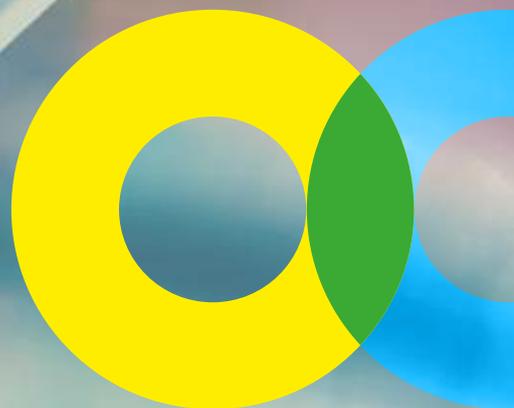

EUROTUBO® slides

Standard size: 26 x 76 mm (± 0.2 mm)

Thickness: 1.1 mm (± 0.1 mm)

Clear, **pre-washed glass**. Ideal for routine analysis.

Displayed in cases of 50 units.

Each case is vacuum flow-packed to protect them from dirtiness.


Each case includes a bag of silica salt. With improved matte band for greater comfort of use: it offers better writing and prevents the slides from sticking together.

Available with polished or cut edges of 45° and 90°.

Expiry date: 24 months.



	code	edges	corners	matte area	case quantity	case weight	case volume
	A D100001	cut 90°	no clipped	no	50 x 50	14.25	0.014
	B D100002	grounded 45°	clipped	no	50 x 50	14.20	0.010
	C D100003	cut 90°	no clipped	yes	50 x 50	14.25	0.014
	D D100004	grounded 45°	clipped	yes	50 x 50	14.25	0.014
NEW	E D100005	grounded 90°	no clipped	no	50 x 50	14.20	0.014
NEW	F D100006	grounded 90°	no clipped	yes	50 x 50	14.20	0.014

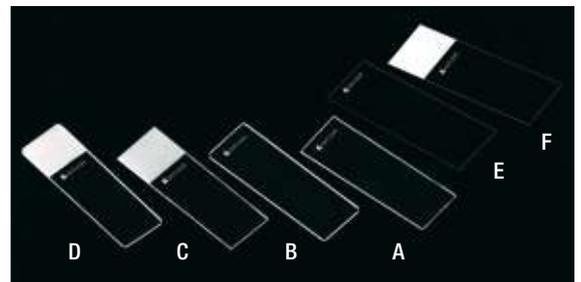

EUROTUBO® cover slides

Standard size: Different according each model. (Please check the table below)

Thickness: 0.145 mm (± 0,015 mm)

Clear, **pre-washed glass**.

Specific use: protection of the sample and lens.

Supplied in cases of a hundred or two hundred units (depending on the code).

Each case includes a bag of silica salt, and is supplied vacuum in tropical pack.

Easy opening, no need for scissors. New improved box which minimize the risk of pack breakage and ensure a better lock.



code	dimensions mm	case quantity	case weight	case volume
D101818	18 x 18	5 x 200	0.30	0.0012
D102020	20 x 20	5 x 200	0.33	0.0011
D102222	22 x 22	5 x 200	0.37	0.0012
D102240	22 x 40	10 x 100	0.65	0.0016
D102424	24 x 24	5 x 200	0.41	0.0011
D102432	24 x 32	10 x 100	0.63	0.0017
D102440	24 x 40	10 x 100	0.68	0.0017
D102450	24 x 50	10 x 100	0.79	0.0017
D102460	24 x 60	10 x 100	0.89	0.0017





STAR FROST® ADHESIVE slides

Standard size: 26 x 76 mm (± 1 mm).

Thickness: 1 mm (± 0,05 mm). Made in Germany.

High definition, **pre-washed glass**.

Hydrophobic. No stain residues, no proteins in the surface.

Marking area: Feature a white, 20 mm. This special printed area is resistant to all common chemicals used in the laboratories, and may be marked with pencil, ballpoint pen, marker and different laboratory-printers (e.g. Leica XPS). Thanks to the single side printing the user always knows the best surface to be used.

Made in Germany.

StarFrost® Adhesive slides do not show striae and disturbing background staining such as the blue or red background obtained with Haematoxylin and Eosin-staining of tissues on Albumin coated glass, or the brown background obtained during the Immunoperoxidase or In-situ-DNS-Hybridisation process. They are ideal for molecular-hybridization, as being DNase and RNase free there is no need of washing them with chromic-acid.

Available with ground edges 90°.

Supplied in boxes of 50 units.

Final package (20 boxes) by **tropical pack**.

Conforms with ISO 8037/1 (european standard). **Expiry date:** 14 months.



Star Frost® Adhesive slides are treated with silane-coating, which ensures cell-and-tissue-adhesion and supports the forming of covalent bonds between the tissue and the glass surface through electrostatic attraction. So tissue and cell substrates stick rapidly on the surface without using further adhesives or protein coating media, saving costs and time. As the underside is untreated, the Star Frost® Adhesive slides do not tend to stick when stacked or transported as with other competitive products.



TROPICAL PACK:

Aluminium film designed to protect slides from humidity.



code	edges	corners	case quantity	case weight	case volume
A100018	grounded 90°	no clipped	20 x 50	5.02	0.004
A100019	grounded 45°	clipped	20 x 50	4.98	0.004

STAR FROST® slides

Standard size: 26 x 76 mm (± 1 mm). Made in Germany.

Thickness: 1 mm.

High definition, **pre-washed glass**.

Hydrophilic (epoxy).

Feature a colour-coded, 20 mm **marking area**.

This special colour-printed area is resistant to all common chemicals used in the laboratories, and may be marked with pencil, ballpoint pen, marker and different laboratory-printers (e.g. Leica XPS). Thanks to the single side printing the user always knows the best surface to be used.



Ground edges at 45° (angle) and clipped corners.

Supplied in boxes of 50 units.

Final package (20 boxes) by tropical pack. **Expiry date:** 14 months.



code	strip colour	case quantity	case weight	case volume
A100010	white	20 x 50	4.94	0.003
A100011	green	20 x 50	4.98	0.004
A100013	blue	20 x 50	4.94	0.004
A100014	yellow	20 x 50	4.94	0.003

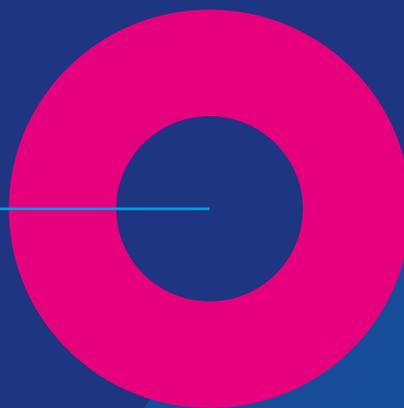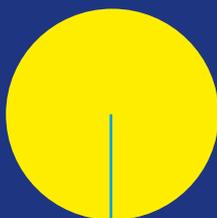

000 «Диаэм»

Москва

ул. Магаданская, д. 7, к. 3 ■ тел./факс: (495) 745-0508 ■ sales@dia-m.ru

www.dia-m.ru

С.-Петербург

+7 (812) 372-6040
spb@dia-m.ru

Новосибирск

+7 (383) 328-0048
nsk@dia-m.ru

Воронеж

+7 (473) 232-4412
vrn@dia-m.ru

Йошкар-Ола

+7 (927) 880-3676
nba@dia-m.ru

Красноярск

+7 (923) 303-0152
krsk@dia-m.ru

Казань

+7 (843) 210-2080
kazan@dia-m.ru

Ростов-на-Дону

+7 (863) 303-5500
rnd@dia-m.ru

Екатеринбург

+7 (912) 658-7606
ekb@dia-m.ru

Кемерово

+7 (923) 158-6753
kemerovo@dia-m.ru

Армения

+7 (094) 01-0173
armenia@dia-m.ru

