


# Automatic Heat Sealer FS-FH10

-Fast, reliable and quiet automated plate sealing solution for any lab





## Features

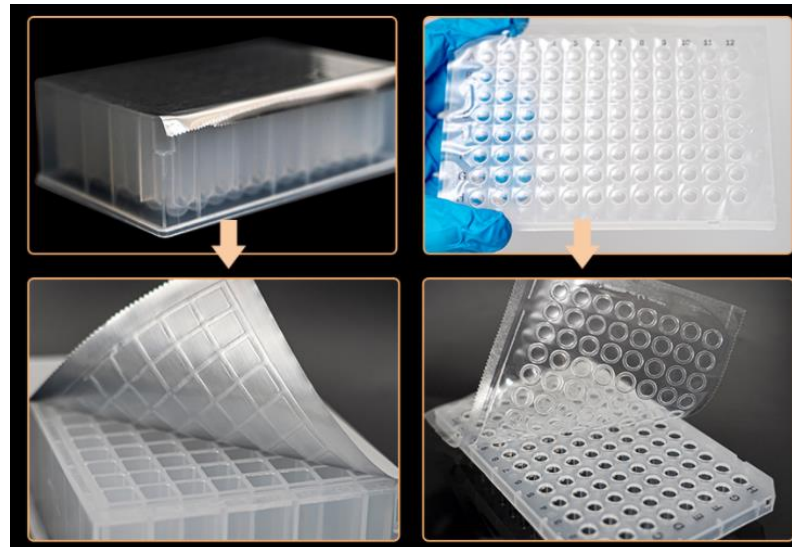1. Fast Sealing: sealing up to **6~8 plates per minute**.
2. Compatible with multiple plates, **adapter can be customized**.
3. Compatible with heatable aluminum foil, and transparent heat films.
4. Tight sealing, no evaporation, good consistency, and repeatability.





Cat No.	Description	Speed	Sealing Temperature Range
FS-FH10	Automatic Heat Sealer	10 seconds per plate	80°C to 200°C

Film	Sealing Conditions	Plate	Photo
Alu foil film (0.1mm thickness)	0 seconds at 180° C to 185° C	2.2ml or 1.6ml Deep well plate	
Transparent film	2 seconds at 145° C to 155° C	PCR 96 well or Elisa plate	




Films: 200 meters per roll can seal more than 1900 plates



**Servicebio**

**Automatic Heat Sealer FS-FH10**

2.2ml Deep well plate : no need adapter

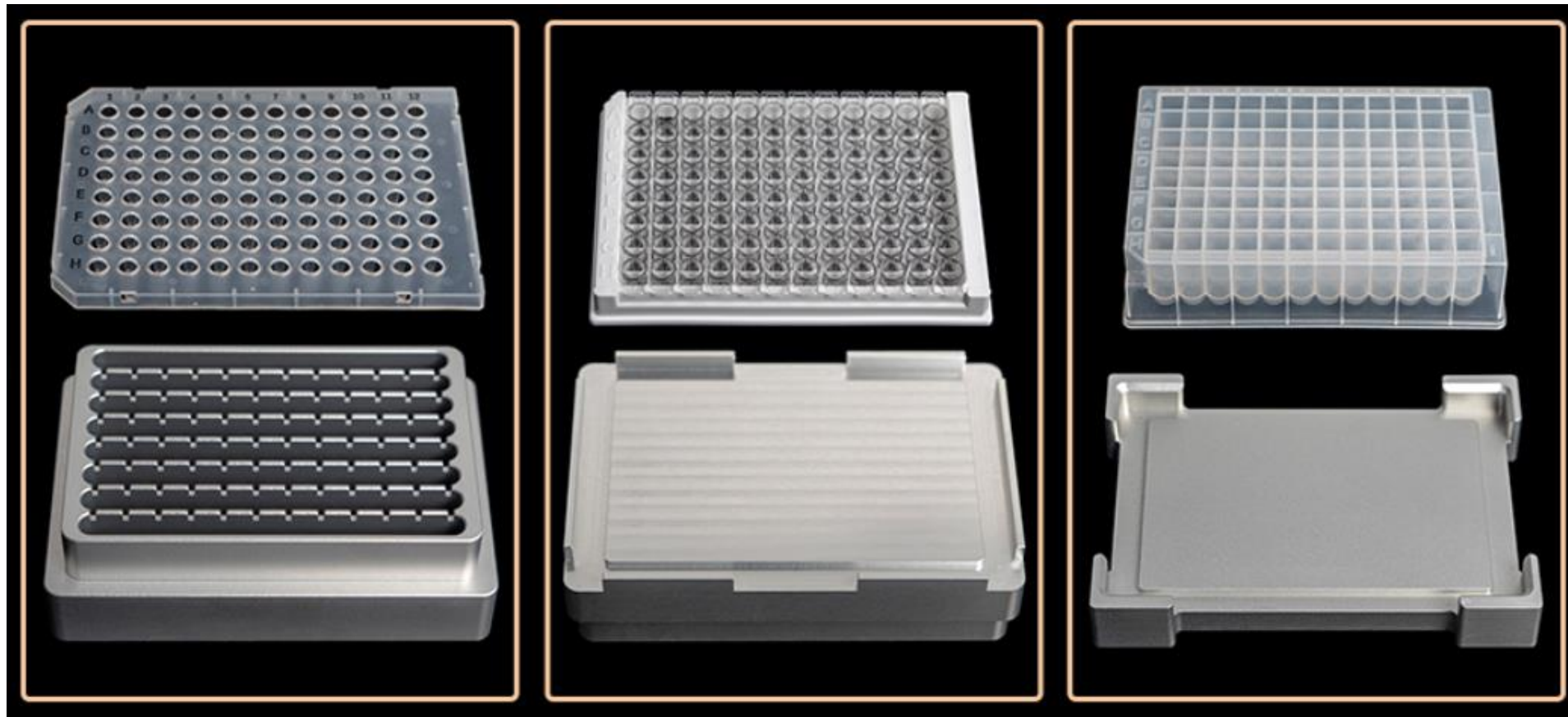

Adapter for 96 well  
PCR plate

Adapter for Elisa plate

Adapter for 1.6ml  
Deep well plate

ЛАБОРАТОРНОЕ ОБОРУДОВАНИЕ ■ ПЛАСТИК ■ СТЕКЛО ■ РЕАКТИВЫ ■ НАБОРЫ

[www.dia-m.ru](http://www.dia-m.ru)  
заказ онлайн